

The Pill by **REVO**TM



Register your Massage Gun for
EXTENDED WARRANTY



revomadic.com/warranty

Massage Gun User Manual

READ AND SAVE THESE INSTRUCTIONS

Inspect Unit Before Use

Contact support@revomadic.com immediately
for any damaged or defective units.

The Pill by REVO™

HAVE AN ISSUE WITH YOUR REVO?

Enjoy hassle-free replacements, If your REVO is defective or damaged in transit, get in touch with us immediately and we'll ship out a new one for **ABSOLUTELY FREE!**

 support@revomadic.com

 **516-450-5170**

 **M-F 9am - 5pm PST**

Scan for
Helpful Resources



To check out more of the products we offer visit:
www.revomadic.com

IMPORTANT TIP

Charge fully before using the device for the first time. **To turn on**, Press & Hold the Power Button. To change Speeds, click the power button until the desired speed is reached. There are 5 speeds.


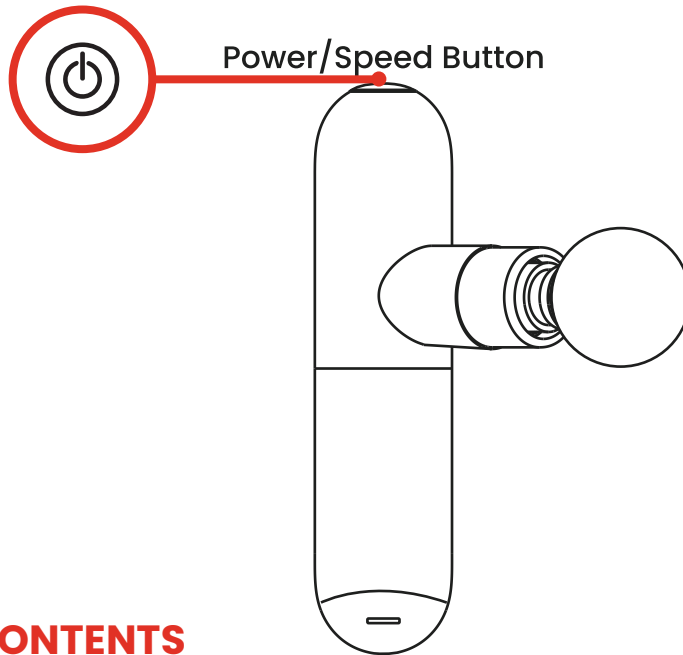



TABLE OF CONTENTS

- | | |
|---------------------|--------------------------|
| 1 Components | 4 Troubleshooting |
| 2 Safety | 5 Specifications |
| 3 Operation | 6 Maintenance |


1. Components

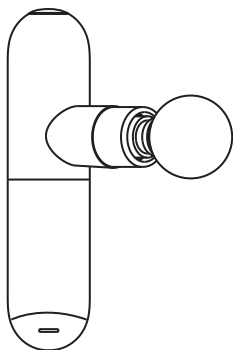


1 Massage Gun

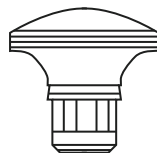


2 Flat (Silicone)

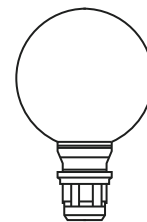


3 50 Ball (EVA foam)

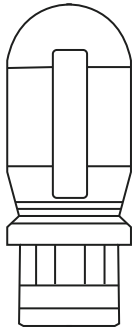

4 Bullet

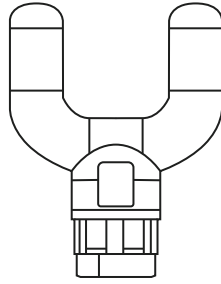

5 Fork

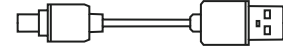

6 Type-C
Charging cable

2. Safety

Caution

Warning:

To reduce the risk of serious injury, read all precautions and instructions in this manual before you use this product. Revomadic assumes no responsibility for personal injury or property damage sustained by or through the use of this product.

- 1 For adults only.** Children or pregnant women are not allowed to use this device. Please ask medical professionals for advice before using this device if carrying diseases or injured. People with physical weakness or mental disorders or slow response shall be under guidance of a guardian. Children should be supervised not to play with this device as a toy.
- 2 When using the device use it only on dry,** clean surfaces of the body by lightly pressing and moving across the body for approximately 60 seconds per region.
- 3 Do not leave it unattended** while charging or operating.
- 4 Do not submerge** the device under water.
- 5 The vents should not be blocked** while operating.
- 6 If the device is heating, stop using it and let it cool down.** If it is overheating, ask a professional to inspect the device before using it again.
- 7 Do not modify the device** by yourself.
- 8 To reduce the risk of electric shock,** do not touch the charge cable or the charging port with wet hands.
- 9 The device should only be used on soft tissue areas of the body** and should not be used on any hard or bony areas of the body, including the head and neck.
- 10 Do not use the device** if you are pregnant, injured, or sick.
- 11 Do not disassemble, crush, burn or short-circuit the battery.**

- 12 When using the device, keep fingers, toes, hair, and other body parts away from the shaft to avoid pinching or entanglement.



Do not put fingers near the head while running.



Do not touch this area to avoid pinching.

Important:


Do not use this product in any of the following cases without the consent of your physician.

Pregnancy, diabetic complications (such as neuropathy or retinal damage), wearing a pacemaker, epilepsy or migraine, disc herniation, advancement of the spine, vertebral detachment, or spinal joint stiffness, recent joint replacement surgery or IUD, having metal needle or steel plate, and any problems related to your health. These contraindications mentioned above do not mean that you are not allowed to use this product, we recommend you to consult your physician before use.

3. Operation

Control Panel

Power/Speed Button

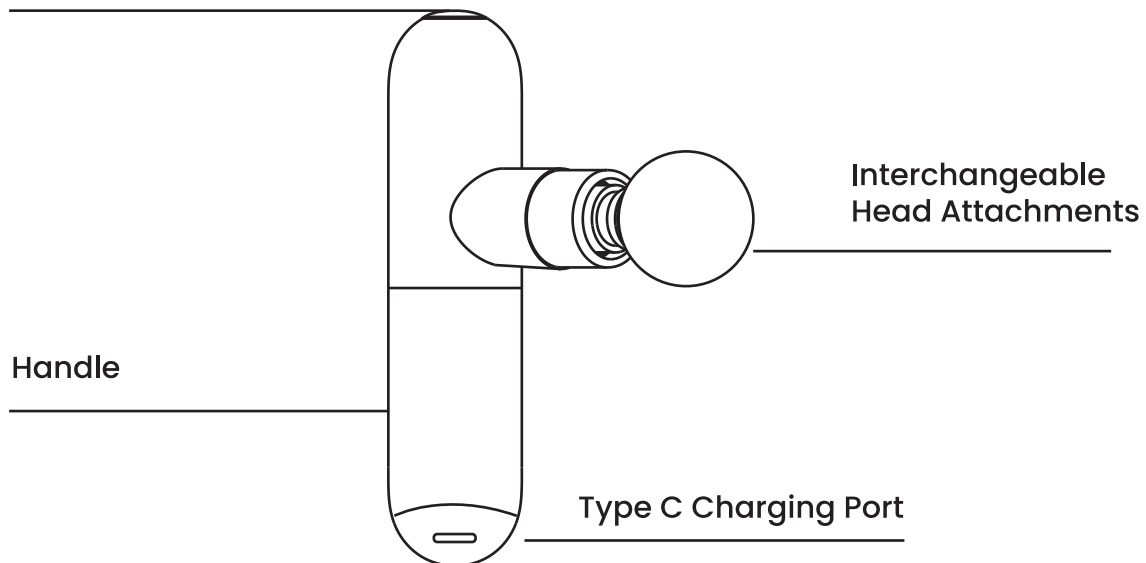

Power On

Hold the Power Button at the Top of the device for 1-2 seconds or until it lights up. Once turned on, the device will be in "Standby" mode.

Click the power button to cycle through the 5 speeds.

Speed Settings

Choose your desired speed (1-5) by pressing the power/speed button.

The device will be in standby after the 5th click.


To protect your bones and fascia, do not use the device for prolonged massage. The device is equipped with 10 minute auto-shutoff feature. Do not use on any area longer than 5 minutes.

Power off

To power off, simply hold the power/speed button down for two seconds.

Battery Level Indicator (LED BAND)

When the battery is fully charged, the battery level indicator will turn green. If the ring is red, then please charge the device.

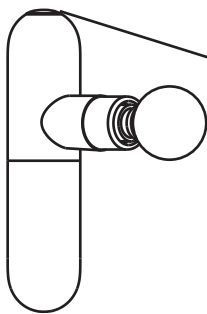


**Power/Speed Button:
Battery Level Indicator (Led Band)**

Using/Replacing Head Attachments

To use a head attachment, insert the desired head directly into the shaft/opening while pushing firmly. To replace, turn the device off then simply rotate or pull out the originally attached head.

Interchangeable Head Attachments

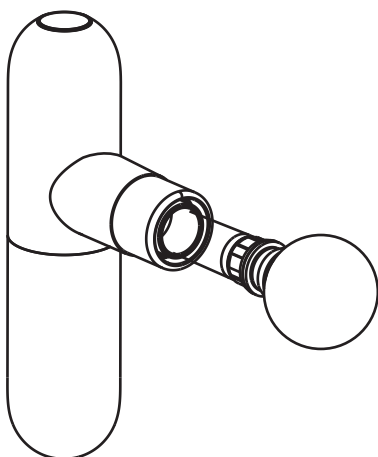

Flat
(silicone)



Φ50 ball
(EVA foam)



Bullet
(Tough PP)



Fork
(Tough PP)

4. Troubleshooting

Troubleshooting


If you experience malfunctions please perform the following checks prior to contacting customer support. Do not disassemble the device or the warranty is voided.

Basic Failure	Possible Causes	Troubleshooting
Cannot be turned on	Dead battery	Charge the device until fully charged
The massaging strength is weak	Battery is low	Charge the device until fully charged.

If all else fails, please get in touch with our support for expert support and replacements.

4. Specifications

Specifications

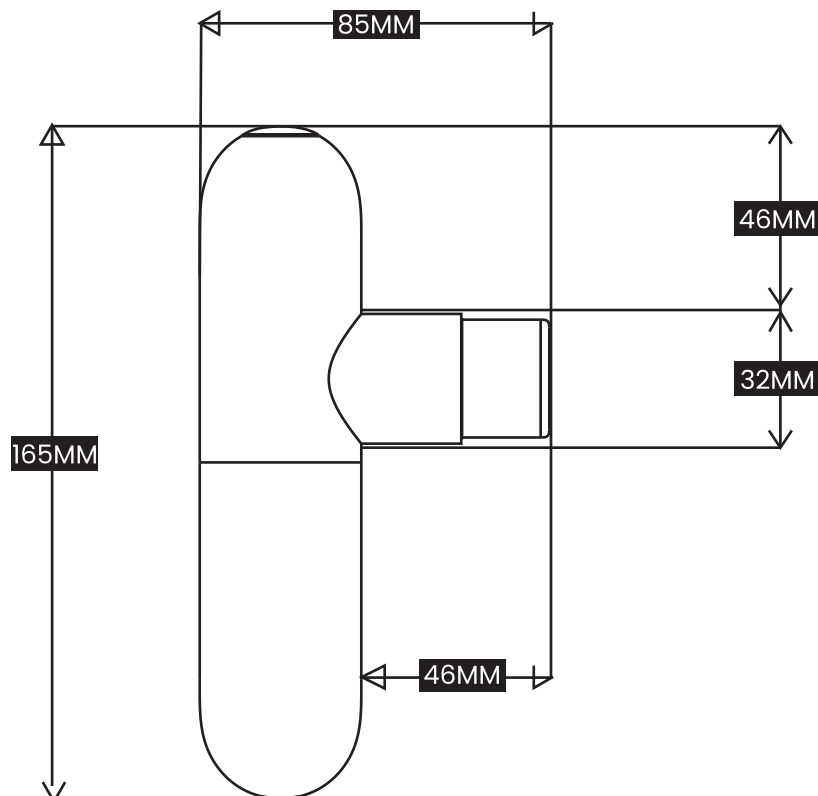

- **Product Size:** _____ 165*85*40mm
- **Product Weight:** _____ 408g
- **Chassis Material:** _____ Aluminum Alloy
- **Amplitude:** _____ 8mm
- **Charging Voltage & Current:** _____ 5V/2A
- **Charging Time:** _____ 3H
- **Stall Force:** _____ $\geq 120\text{N}$
- **Rotation Speed Levels:** _____ 2000/2300/2600/2900/3200rpm
- **Cell Capacity:** _____ 1300mAh
- **Battery Specification:** _____ Lithium Polymer Battery (non-explosive)

6. Maintenance

Storage and Maintenance

Storage

If the device has been stored for more than 3 months, inspect and recharge it before use. Do not place it on the floor, next to walls, or on the roof. Make sure ventilation is good and keep the device away from direct sunlight and corrosive gases.

Transportation

The Pill is safe and permitted for international air travel in checked in or carry on suitcases or backpacks.

Power Supply Check

Check the power supply before charging the device. It should not exceed 5V/2A.

Maintenance

This device should be protected from moisture, rodents, and pollution. Please turn off the device when cleaning. After use, please put the device in a safe place. Do not place heavy objects on the unit. When the device is malfunctioning, please ask professional technicians to repair it, do not disassemble it by yourself.

Disposal & Recycling

This device should not be disposed of together with household garbage. When you decide to dispose of this device, please follow the local environmental laws and disposal guidelines.